

V61

He-Ne Laser

Fritz Ali Agildere
fritz.agildere@udo.edu

Jan Lucca Viola
janlucca.viola@udo.edu

Experiment: November 25, 2024

Submission: November ??, 2024

TU Dortmund – Department of Physics

Contents

1	Objective	1
2	Background	1
2.1	Components of a laser	1
2.2	Processes in the active medium	2
2.3	Necessity of multiple levels	4
2.4	Stability for different resonators	6
2.5	Transverse and longitudinal modes	8
2.6	Doppler broadening of the transition	9
2.7	Brewster window and polarization	9
2.8	Wavelength from the diffraction pattern	10
3	Procedure	10
3.1	Aligning the Laser	10
3.2	Verifying the Stability Condition	10
3.3	Observing Transverse Modes	10
3.4	Determining the Polarization	11
3.5	Analyzing Spectra in Multimode Operation	11
3.6	Measuring the Wavelength	11
4	Results	11
5	Discussion	11
	References	11
	Appendix	12

1 Objective

To understand the characteristics of a He-Ne laser, several different configurations are adjusted and their radiation properties measured. This includes the wavelength, intensity distribution, polarization, mode spectrum, as well as the influence of mirror type and resonator length.

2 Background [1]

Lasers are ubiquitous tools of modern physics due to their useful properties, characterized by the emission of coherent light with narrow spectral linewidth, low divergence and high power density. They are named after the acronym for light amplification by stimulated emission of radiation, describing the fundamental mechanism for the production of laser radiation. This will be explored in the following, both in the general as well as the special case of the He-Ne laser.

2.1 Components of a laser

The basic setup of a typical laser consists of three main components, namely an active medium, a pumping mechanism and the resonator cavity.

Inside the active medium, realized using materials such as semiconductors or gas mixtures, photons are emitted from atomic transitions to energetically lower states. The energy difference ΔE between the involved electron levels is therefore the main determinant of wavelength λ and frequency f via $\Delta E = hf$.

To excite electrons in the active medium to higher levels, an energy source is required. This is the role of the pumping mechanism, which can be implemented using electrons or photons. The latter case is called optical pumping, as another separate light source tuned to the respective ΔE value is used to induce transitions to excited states.

Amplification of the emitted radiation is achieved in the active medium. Instead of using superradiant lasers, which have high gain factors and divergence, or impractically long constructions, mirrors can be used to create a resonator cavity. The resulting standing waves correspond to multiple passes through the material and can generate a stable beam with low divergence, which can exit through a semitransparent window. The mirror geometry can be adapted to the desired function with flat or concave designs.

2.2 Processes in the active medium

There are three main processes shown in Figure 1 occurring inside the active medium to facilitate the operation of a laser.



Figure 1: Schematic depiction of relevant processes inside the active medium. [1]

Raising the energy of an electron by $\Delta E = E_2 - E_1$ requires the annihilation of an incident photon that fulfills the condition

$$\Delta E = hf_{12} ,$$

where h is the Planck constant. This process is referred to as absorption. The number of transitions per time and volume is proportional to the density of ground state electrons N_1 as well as the photon flux or number per area and time φ via

$$\left. \frac{dN_1}{dt} \right|_{\text{ab}} = -\sigma_{12} N_1 \varphi ,$$

with σ_{12} denoting the effective cross section for absorbing a photon. From this also follows the typical exponential intensity reduction

$$\left. \frac{dI}{dx} \right|_{\text{ab}} = -\sigma_{12} N_1 I ,$$

where $\alpha = -\sigma_{12} N_1$ gives the absorption coefficient.

When an atom is in an excited state, it returns to the ground state after a time interval, the duration of which follows some random distribution with mean lifetime τ . Due to its stochastic nature, this process is called spontaneous emission and has no predefined direction or phase. The density in the higher level then follows

$$\left. \frac{dN_2}{dt} \right|_{\text{sp}} = -\tau^{-1} N_2 .$$

Besides this, emission can also be initiated by an incoming photon of appropriate frequency. This is called stimulated emission and results in the production of radiation with the same energy, direction and phase as the inducing quantum. As the inverse process to absorption,

$$\left. \frac{dN_2}{dt} \right|_{\text{st}} = -\sigma_{21} N_2 \varphi ,$$

describes the time derivative and

$$\left. \frac{dI}{dx} \right|_{\text{st}} = \sigma_{21} N_2 I ,$$

the corresponding intensity relation. This means that stimulated emission leads to an increase in intensity, serving as a potential mechanism for amplification when there are more electrons in the excited state than in the ground state and losses are compensated for. This phenomenon is referred to as population inversion.

The cross sections can be identified with the Einstein coefficients B_{ij} via

$$\sigma_{ij} = B_{ij} h f_{ij} / c ,$$

where c is the speed of light in vacuo. Furthermore, for emission and absorption, the thermodynamic or quantum mechanical relation

$$g_1 \sigma_{12} = g_2 \sigma_{21}$$

holds, with g_1 and g_2 defining the degrees of degeneracy for the ground and excited states. Hereafter, it is assumed that E_1 and E_2 have the same number of sublevels, so $g_1 = g_2$ for $\sigma_{12} = \sigma_{21}$ and

$$B \equiv B_{12} = B_{21} .$$

The reciprocal decay timescale defines another Einstein coefficient

$$A = \tau^{-1} ,$$

with which the stationary spectral radiance

$$\rho_s \equiv \frac{A}{B} = \frac{8\pi h f_{12}^3}{c^3}$$

can be written. Introducing the general spectral radiance

$$\rho = \varphi h f_{12} / c$$

and requiring $N = N_1 + N_2$ to be constant for a system of two energy levels, one finds

$$\frac{dN_1}{dt} = \left. \frac{dN_1}{dt} \right|_{\text{ab}} - \left. \frac{dN_2}{dt} \right|_{\text{st}} - \left. \frac{dN_2}{dt} \right|_{\text{sp}} = \rho B (N_2 - N_1) + A N_2 = -\frac{dN_2}{dt} .$$

For $\Delta N = N_2 - N_1$ then follows that

$$\frac{d\Delta N}{dt} = -2\frac{dN_1}{dt} = -2\rho B\Delta N - 2AN_2 + AN_1 - AN_1 = -2\rho B\Delta N - A\Delta N - AN.$$

After some time an equilibrium is reached inside the active medium, resulting in a vanishing time derivative. In this case, solving for the stationary number difference

$$\Delta N_s = -\frac{AN}{A + 2\rho B} = -\frac{N}{1 + 2\rho/\rho_s}$$

yields $\Delta N_s < 0$ for any system with only two energy levels.

2.3 Necessity of multiple levels

This result directly contradicts the requirement of population inversion $\Delta N_s > 0$ necessary for the amplification through stimulated emission as discussed previously, preventing the usage of two level systems as the active medium. Adding more energy levels as depicted in Figure 2 solves this problem. Instead of immediately relaxing back to the ground state via spontaneous emission, excited electrons now decay very quickly from E_3 to E_2 and E_1 to E_0 while the E_2 to E_1 transition takes longer. This means that $A_{21} < A_{32}$ as well as $A_{21} < A_{10}$ and results in a distribution similar to what is shown below.



Figure 2: Exemplary energies and population densities for multiple levels. [1]

One then expects $N_0, N_2 \gg N_1, N_3$ for $N \approx N_0 + N_2$ and $\Delta N \approx N_2$ in the stationary configuration. Accordingly, a population inversion $\Delta N_s > 0$ is trivial to achieve, making four level systems a suitable choice for laser construction.

Such a system is realized by a He-Ne laser when the red mode at $\lambda = 633\text{ nm}$ is used. Table 1 indicates that this corresponds to a transition from the $3s_2$ to the $2p_4$ level, on which Figure 3 provides more detailed information.

Table 1: Properties for different transitions of the He-Ne laser. [1]

Color	Wavelength (nm)	Transition (Paschen notation)	Output power (mW)	Spectral width (MHz)	Gain (%/m)
Infrared	3391	$3s_2 \rightarrow 3p_4$	>10	280	10,000
Infrared	1523	$2s_2 \rightarrow 2p_1$	1	625	
Infrared	1153	$2s_2 \rightarrow 2p_4$	1	825	
Red	640	$3s_2 \rightarrow 2p_2$			
Red	635	$3s_2 \rightarrow 2p_3$			
Red	633	$3s_2 \rightarrow 2p_4$	>10	1500	10
Red	629	$3s_2 \rightarrow 2p_5$			
Orange	612	$3s_2 \rightarrow 2p_6$	1	1550	1.7
Orange	604	$3s_2 \rightarrow 2p_7$			
Yellow	594	$3s_2 \rightarrow 2p_8$	1	1600	0.5
Green	543	$3s_2 \rightarrow 2p_{10}$	1	1750	0.5

In this setup, the active medium is a mixture of helium and neon gases. Helium atoms are excited to metastable states via electric discharge before colliding with neon atoms to provide excitation and transfer kinetic energy via



The excess energy $E \simeq 100$ meV is dissipated as heat after the resonant transfer, as it measures about two times the 300 K thermal energy. Due to the selection rules, the upper 2s and 3s levels have lifetimes of the order 100 ns because they can only decay to p levels, while lower states exhibit shorter 10 ns timescales.



Figure 3: Energy level diagram of a He-Ne laser in Paschen notation. [1]

2.4 Stability for different resonators

To characterize spherical mirrors, the parameter

$$g_k = 1 - \frac{L}{R_k}$$

presents a convenient alternative to the curvature radii. With this, the radius w_0 at the waist of a Gaussian beam as depicted in Figure 4 can be calculated via

$$w_0^4 = \frac{\lambda^2 L^2 g_1 g_2 (1 - g_1 g_2)}{\pi^2 (g_1 + g_2 - 2g_1 g_2)^2}.$$

It is immediately obvious that for w_0 to be a real solution,

$$0 \leq g_1 g_2 \leq 1$$

needs to hold, defining a stability condition for any mode to exist in a given resonator.

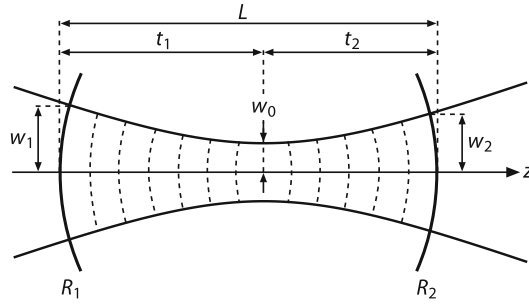


Figure 4: Adaptation of a Gaussian beam to a spherical mirror resonator. [1]

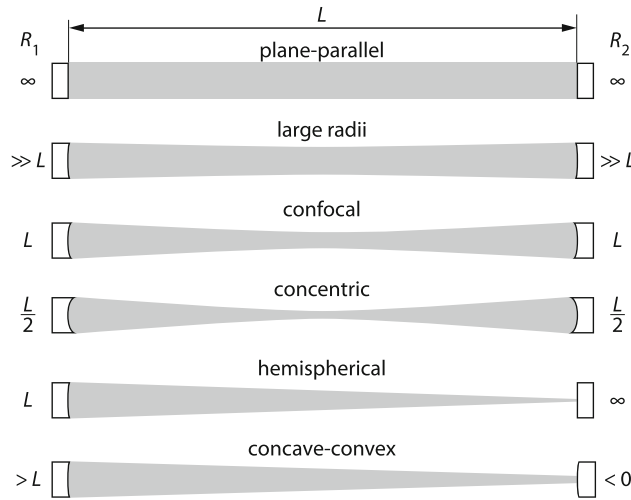


Figure 5: Types of stable resonator configurations. [1]

Table 2: Available mirrors for the construction of the cavity.

Mirror	Design	Surface
planar	$R_1 = \infty$ flat / flat	HR (high reflectivity) $R \geq 99\%$
concave	$R_1 = 1000$ mm spherical / flat	HR (high reflectivity) $R \geq 99\%$
concave	$R_1 = 1400$ mm spherical / flat	HR (high reflectivity) $R \geq 99\%$
concave	$R_2 = 1400$ mm spherical / flat	OC (out coupling) $T \simeq 2\%$

Figure 5 displays some realizable mirror configurations, while Table 2 lists available parts for the setup at hand. Their stability factors are evaluated graphically in Figure 6 to determine the maximum possible resonator length.

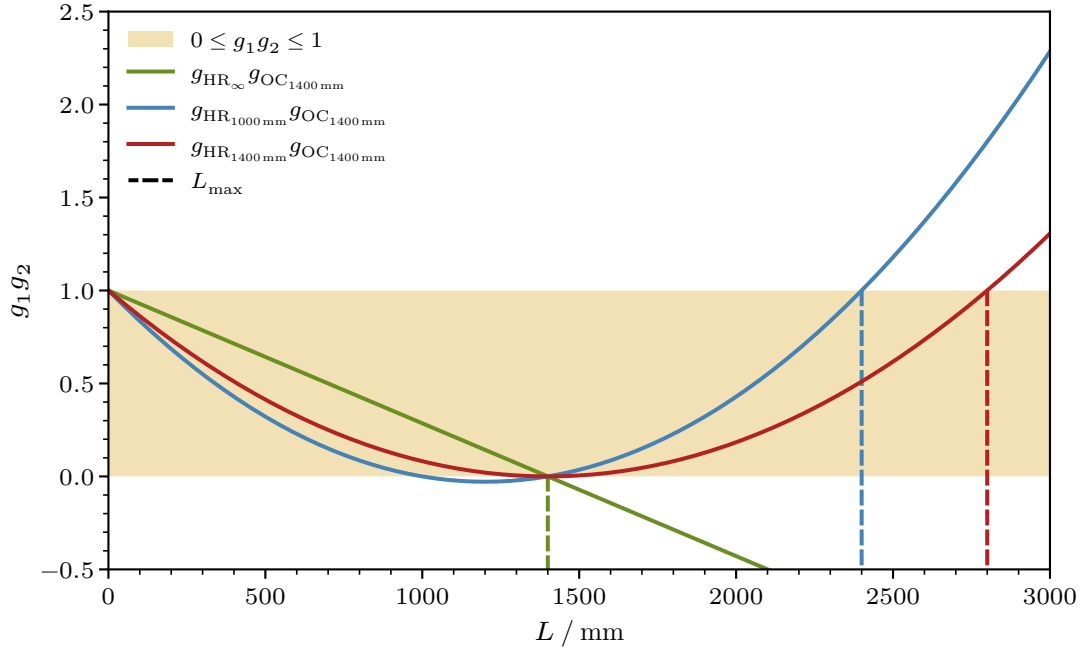


Figure 6: Stability diagram as functions of the resonator length.

It should be noted that an unstable resonator with $g_1g_2 < 0$ or $g_1g_2 > 1$ does not necessarily preclude laser operation. Unstable lasing implies large diffraction losses, as significant parts of laser power leak around the mirror edges, in turn requiring high optical gain in the active medium to compensate. The resulting configuration can however be less sensitive to disturbances such as thermal lensing or misalignment than a stable resonator, providing desirable properties depending on the application.

2.5 Transverse and longitudinal modes

During normal operation, different transverse electromagnetic modes (TEM) interfere in the plane perpendicular to the beam axis. Their orders are denoted using indices like $\text{TEM}_{r\vartheta}$ and TEM_{xy} to give the integer number of nodes in direction of the respective coordinate. Assuming a Gaussian beam profile, circular symmetry can be parametrized by multiplication with Laguerre polynomials. If modes are selected using a thin wire as is the case here, a rectangular symmetry is required. A cut of the electric field strength in x direction then follows

$$E_n(\xi) = H_n(\xi)e^{-\xi^2/2},$$

where $H_m(\xi)$ is the appropriate Hermite polynomial and $\xi = \sqrt{2}x/w_0$ connects the parameter to the spacial coordinate. The first three solutions

$$H_0 = 1, \quad H_1(\xi) = 2\xi, \quad H_2(\xi) = 4\xi^2 - 2$$

lead to the patterns shown in Figure 7 below.

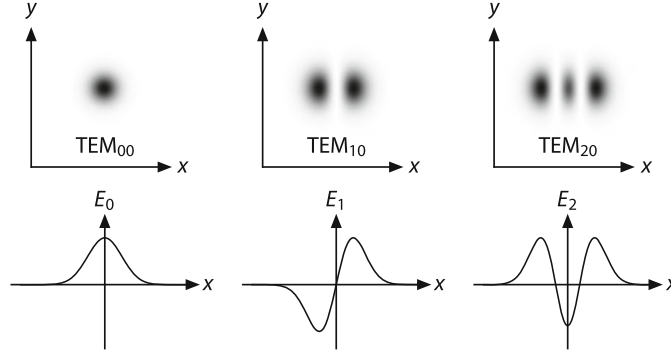


Figure 7: Observable TEM_{xy} modes for rectangular geometry. [1]

Obtaining the intensity, which is the quantity that is actually measured, requires squaring the absolute value of the electric field amplitude. From $I \propto |E|^2$ follows

$$\begin{aligned} I_{\text{TEM}_{00}}(x) &= I_{00}e^{-2x^2/w_0^2}, \\ I_{\text{TEM}_{10}}(x) &= I_{10}x^2e^{-2x^2/w_0^2}/w_0^2, \\ I_{\text{TEM}_{20}}(x) &= I_{20}(4x^2/w_0^2 - 1)^2e^{-2x^2/w_0^2} \end{aligned}$$

as the theoretically expected distributions with arbitrary normalization factors. Additional degrees of freedom may be introduced to account for distortive effects such as a finite background or an asymmetric shape due to inaccurate centering of the probe.

Longitudinal modes in direction of beam expansion instead interfere to produce beating patterns in the signal. This is due to slight frequency variations that result from the multiple allowed standing wave solutions inside the resonator.

2.6 Doppler broadening of the transition

Lorentzian

Gaussian

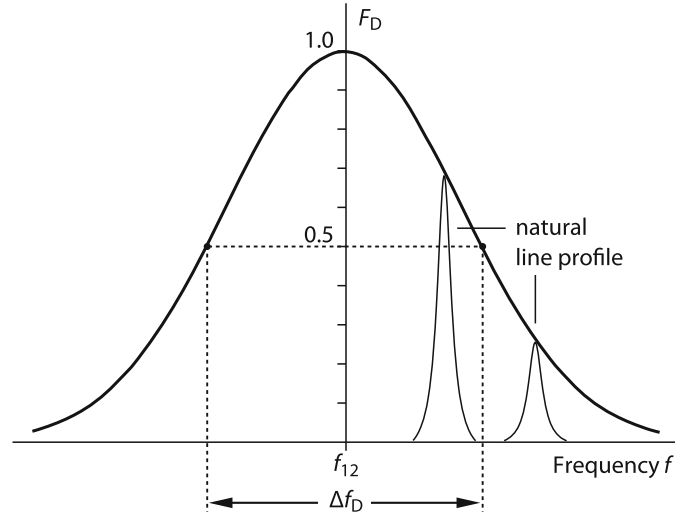


Figure 8: Typical Doppler broadening compared to natural line width. [1]

2.7 Brewster window and polarization

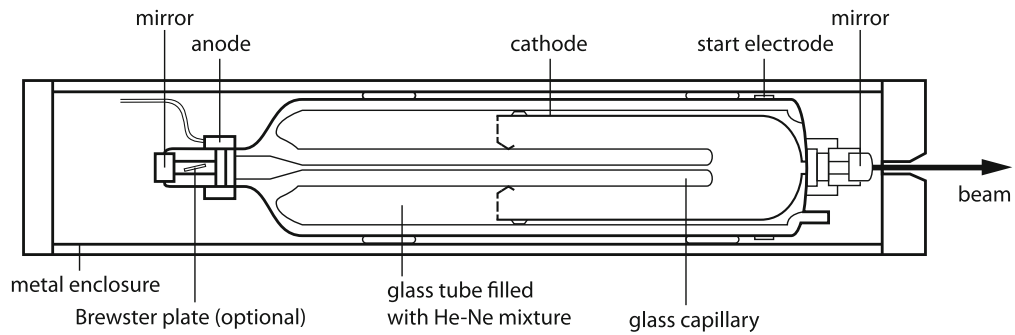


Figure 9: Schematic setup of a He-Ne laser generating polarized output. [1]

2.8 Wavelength from the diffraction pattern

3 Procedure

3.1 Aligning the Laser

The alignment of the He-Ne laser is initiated by positioning an auxiliary alignment laser (wavelength: 532 nm, maximum power: 1 mW, reduced power: 0.2 mW) on the optical bench. A target screen with a cross mark is placed at the end of the optical rail. The alignment laser should be adjusted such that its beam passes through the center of the cross on the target screen. This ensures that the alignment laser coincides with the optical axis of the He-Ne laser.

Next, the He-Ne laser components are positioned in the following order: the laser tube (length: 408 mm, diameter: 1.1 mm), resonator mirrors (diameter: 12.7 mm), and Brewster windows. These components together form the laser resonator, with the Brewster windows ensuring minimal loss while defining the polarization direction. The alignment of these components is critical: the alignment laser's back reflections should be made to hit the target screen's cross at the center, indicating proper alignment along the optical axis.

3.2 Verifying the Stability Condition

Once the laser has been aligned, the stability condition must be verified. The laser is adjusted to its maximum power using a photodiode, and the maximum resonator length is set by gradually increasing the gap between the two resonator mirrors. Throughout the process, the laser power is continuously readjusted. With a well-aligned setup, the system should approach the theoretical value from the stability condition. This step is repeated for different resonator lengths to study the effect on the laser's stability.

3.3 Observing Transverse Modes

To observe transverse electromagnetic (TEM) modes, a thin tungsten wire (diameter: 0.005 mm) is placed between the resonator mirror and the laser tube. This wire stabilizes different modes, which can be observed on an optical screen. A scattering lens may be used to enlarge the laser beam, making it easier to identify the modes. The wire functions to stabilize the laser beam, enabling clearer mode identification. The optical screen is then replaced with a photodiode to measure the intensity distribution for at least two modes. The measured intensity distributions are plotted and compared with theoretical expectations to validate the mode stability.

3.4 Determining the Polarization

The laser's polarization is determined by placing a polarizer behind the outcoupling mirror. The intensity of the laser beam is measured with a photodiode as the polarizer is rotated. The Brewster windows, which minimize reflection losses, ensure a well-defined polarization direction. By comparing the experimental intensity distribution with theoretical calculations, the polarization characteristics of the laser are evaluated, and the influence of the Brewster windows and resonator mirrors on the polarization is examined.

3.5 Analyzing Spectra in Multimode Operation

In the absence of a Fabry-Perot etalon, the laser operates in multimode, meaning several longitudinal modes coexist. This leads to temporal intensity variations due to the beating between modes. To analyze these beat frequencies, a fast photodiode with a bandwidth up to 1 GHz is used, and the Fourier spectra are recorded for various resonator lengths with a spectrum analyzer. The spread of the neon transition is compared with the distance between the longitudinal modes, and the multimode operation is justified. Additionally, the dependence of the beat frequency on the resonator length is investigated.

3.6 Measuring the Wavelength

The wavelength of the He-Ne laser is determined by using diffraction patterns produced by a slit and diffraction grating. The diffraction maxima and minima are measured to accurately determine the wavelength. This method provides a precise measurement of the laser's wavelength and can be used to verify the laser's output characteristics.

4 Results

5 Discussion

References

- [1] Hans Joachim Eichler, Jürgen Eichler, and Oliver Lux. *Lasers. Basics, Advances and Applications*. Springer Cham, 2018. ISBN: 978-3-319-99895-4. DOI: <https://doi.org/10.1007/978-3-319-99895-4>.

Appendix